

COVID-19 Test Home Collection Kit



Let's get started.

WWW.PIXEL.LABCORP.COM

Kit includes:




1 Saline
collection tube



1 foam swab



1 Biohazard
specimen bag

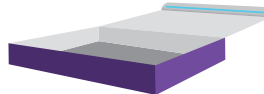

1 Pixel box



FedEx return pack
UN 3373

Intended use: to collect a nasal swab for purpose of in vitro diagnostic testing in a CLIA laboratory.
Rx only

IVD

Distributed by Pixel by LabCorp
531 South Spring Street
Burlington, NC 27215

This home collection kit has not been FDA cleared or approved.

This home collection kit has been authorized by FDA under an EUA.

This home collection kit has been authorized only for the home collection and maintenance of nasal swab specimens as an aid in detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens.

This home collection kit in combination with the authorized test is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostic tests for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Act, 21 U.S.C. § 360bbb-3(b)(1), unless the authorization is terminated or revoked sooner

PXC62920

Self-Collection kit instructions

For In Vitro Diagnostic Use



1

Please read all of the instructions before beginning collection. You may also watch our video on www.pixel.labcorp.com to view an example of a successful collection.

REGISTER ONLINE

2

When you are ready to collect your sample, register your kit online at www.pixel.labcorp.com/register and type in the 12-digit registration code that is located on your collection tube, Pixel box, and registration card.




3

Visit <https://bit.ly/fedexdropbox> to view FedEx drop box locations and pickup schedules. It's important to bring your sample to a drop box on the same day you collect it before the last Express pick up. Or refer to the flyer in your kit to schedule a FedEx pickup on the day of collection.




Wash and dry hands before opening the kit. Open your kit and place all the contents on a clean, dry surface.

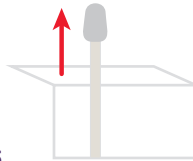


5

Take the foam swab out of its package. Do not touch the tip of the swab with your hands.

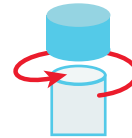


6

Screw off the top of the collection tube. Hold swab in one hand and collection tube in the other being careful not to spill the liquid. Do not drink the liquid.

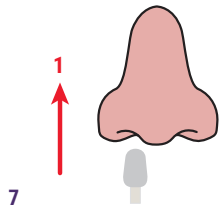


7

Insert the tip of the foam swab into one nostril. The foam swab does not need to be inserted far – insert just until the foam tip of the swab is no longer visible. Rotate the swab in a circle around the entire inside edge of your nostril at least 3 times.

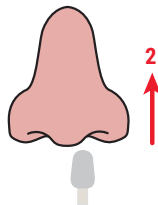

8

Take the foam swab out of your nostril. Using the same end of the foam swab, repeat step 7 in your other nostril.

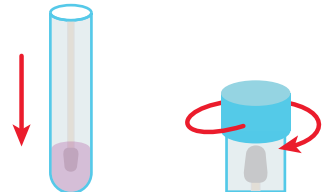


9

Remove the foam swab from your second nostril and place in the collection tube. The end of the foam swab that went into your nose should be placed into the tube first so that it sits down in the liquid. Screw the top of the collection tube back on.



10

Wash and dry hands thoroughly again. Insert collection tube into the biohazard specimen bag. Seal the biohazard specimen bag by closing the zip lock seal.

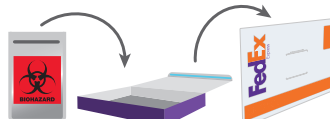

11

Place the specimen bag back into your purple Pixel shipping box and close the lid. Place the shipping box into the FedEx return pack. Remove the adhesive cover strip and seal the return pack closed.



12

Deliver the postage paid, pre-addressed FedEx return pack to a FedEx drop box. It's important to bring your sample to a drop box on the same day you collect it before last Express pick up. Or refer to the flyer in your kit to schedule a FedEx pickup on the day of collection.