

FTD™ SARS-CoV-2

Current Revision and Date 11416299_en Rev. B, 2020-07

Product Name FTD SARS-CoV-2 (FTD-114-96)

REF 11416302



Specimen Types Upper respiratory specimens (such as nasal swabs, nasopharyngeal swabs, oropharyngeal swabs, nasopharyngeal wash/aspirate or nasal aspirate) and Bronchoalveolar lavage (BAL)

Processed Sample Volume 200 µL / 500 µL required (see details in the respective extraction sections)

FTD SARS-CoV-2 was validated with the NucliSENS® easyMAG® (bioMérieux)/Applied Biosystems® 7500 Real-Time PCR System (ThermoFisher Scientific) and the VERSANT® kPCR Molecular System (Siemens Healthineers).

Table of Contents

Intended Use	2
Summary and Explanation.....	2
Principles of the Procedure	3
Method	3
Reagents.....	4
Warnings and Precautions.....	4
Handling Requirements.....	5
Storage and Handling	5
Specimen Collection and Handling	6
Collecting the Specimen	6
Storing and Transporting the Specimen.....	7
Procedure	8
Materials Provided	8
Materials Required but Not Provided	9
Assay Procedures	11
Quality Control	34
Criteria for a Valid Run	34
Results.....	35
Interpretation of Results.....	35
Limitations.....	37
Conditions of Authorization for the Laboratory	38
Performance Characteristics	39
Limit of Detection – Analytical Sensitivity	39
Inclusivity – Analytical Sensitivity.....	41
Cross-Reactivity – Analytical Specificity.....	41
Clinical Performance.....	46
NucliSENS® easyMAG® / Applied Biosystems® 7500.....	46
VERSANT® kPCR Molecular System	47
Troubleshooting	48
Technical Assistance	49
References.....	49
Definition of Symbols.....	50
Legal Information	51

Intended Use

FTD SARS-CoV-2 is a real-time polymerase chain reaction (RT-PCR) test intended for the qualitative detection of nucleic acid from SARS-CoV-2 in upper respiratory specimens (such as nasal swabs, nasopharyngeal swabs, oropharyngeal swabs, nasopharyngeal wash/aspirate or nasal aspirate) and bronchoalveolar lavage from individuals suspected of Coronavirus Disease 2019 (COVID-19) by their healthcare provider. Testing is limited to laboratories – certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, to perform high complexity tests. Results are for the identification of SARS-CoV-2 RNA. The SARS-CoV-2 RNA is generally detectable in upper respiratory specimens and bronchoalveolar lavage during the acute phase of infection. Positive results are indicative of the presence of the SARS-CoV-2 RNA; clinical correlation with patient history and other diagnostic information is necessary to determine the patient infection status. Positive results do not rule out bacterial infection or co-infection with other viruses. Laboratories within the United States and its territories are required to report all positive results to the appropriate public health authorities.

Negative results do not preclude SARS-CoV-2 infection and should not be used as the sole basis for patient management decisions. Negative results must be combined with clinical observations, patient history, and epidemiological information.

FTD SARS-CoV-2 is intended for use by trained qualified laboratory personnel specifically instructed and trained in the techniques of RNA extractions and real-time PCR (*in vitro* diagnostic procedures). FTD SARS-CoV-2 is only for use under the Food and Drug Administration's (FDA) Emergency Use Authorization.

Summary and Explanation

On December 31st, 2019, the World Health Organization (WHO) was informed of multiple cases of pneumonia of unknown etiology detected in Wuhan City, Hubei Province of China. Soon, a new strain of coronavirus, SARS-CoV-2, observed for the first time in humans was identified to be the cause of this new disease later called COVID-19. On January 30th, 2020, WHO declared SARS-CoV-2 as a Public Health Emergency of International Concern. Since its emergence it has rapidly spread worldwide, causing a massive global outbreak, which has reached the status of a pandemic.

The first symptoms of the COVID-19 are not very specific. People may experience runny nose, headache, muscle pain and tiredness. Fever, cough and respiratory signs often occur 2 or 3 days later and can lead to severe pneumonia and death. The risk of developing more severe symptoms of COVID-19 are currently unknown; however, individuals with pre-existing conditions may be more at risk of developing severe symptoms. The duration of incubation is on average 5 days, with extremes of 2 to 12 days.¹

FTD SARS-CoV-2 is an aid in the identification of COVID-19 disease by the detection of SARS-CoV-2 RNA extracted from upper respiratory specimens (such as nasal swabs, nasopharyngeal swabs, oropharyngeal swabs, nasopharyngeal wash/aspirate or nasal aspirate) and bronchoalveolar lavage (BAL) from individuals suspected of COVID-19 by their healthcare provider.

Principles of the Procedure

Method

FTD SARS-CoV-2 is a real-time reverse transcription polymerase chain reaction (rRT-PCR) test for the detection of RNA from SARS-CoV-2 in human upper respiratory and BAL specimens.

Nucleic acids should be first extracted from the specimen types listed in the *Intended Use* section, with addition of the internal control (IC).

The eluate with purified nucleic acids of SARS-CoV-2 is added to a master mix to enable the RT-PCR reaction using Applied Biosystems® 7500 Real-Time PCR Thermocycler (nucleic acids extracted with NucliSENS® easyMAG®) or the VERSANT® kPCR Molecular System. The master mix contains enzyme, buffer, deoxyribonucleotide triphosphate (dNTPs) and synthetic primers and probes specific for the targeted sequences. The advantage of using multiple primer/probe pairs in a single reaction mixture is to simultaneously detect different targets in one reaction.

In the presence of the target, primers and probes will hybridize to the specific sequence and allow amplification by the polymerase. The different probes include fluorescent dyes and quenchers in close proximity to each other, limiting the fluorescence emitted. However, during amplification, the polymerase extends the new strand and degrades the fluorescent dye-labeled probe using its exonuclease activity. This results in the separation of the fluorescent dye from the quencher, thereby allowing the emission of fluorescence.

The level of fluorescence increases with the generation of amplicons and is proportional to the amount of pathogen nucleic acid contained in the sample. The increased fluorescence level is reported as a cycle threshold (Ct) value by the real-time thermocycler.

The assay uses primer and probe sets that target N gene and ORF1ab region of SARS-CoV-2. The mix further includes a primer and probe set to detect a sequence in the genome of equine arteritis virus (EAV) that serves as an internal control (IC).

Reagents

Warnings and Precautions

CAUTION

Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.

- For use under an Emergency Use Authorization (EUA) only.
 - For *in vitro* diagnostic use (IVD).
 - For prescription use only.
 - FTD SARS-CoV-2 has not been FDA cleared or approved; the test has been authorized by FDA under an Emergency Use Authorization (EUA) for use by laboratories certified under the Clinical Laboratory Improvement Amendments (CLIA) of 1988, 42 U.S.C. §263a, to perform high complexity tests.
 - FTD SARS-CoV-2 has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens.
 - FTD SARS-CoV-2 is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of *in vitro* diagnostic tests for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Act, 21 U.S.C. §360bbb-3(b)(1), unless the authorization is terminated or revoked sooner.
 - Laboratories within the United States and its territories are required to report all positive results to the appropriate public health authorities.
-

Safety data sheets (SDS) are available at www.siemens-healthineers.com/sds. Strict adherence to the following warnings and precautions are required when running FTD SARS-CoV-2.



WARNING

The IC contains lysis buffer.



WARNING

FTD SARS-CoV-2 IC:

Hazardous ingredient: Maleic acid (0.1% [w/w])

H317: May cause an allergic skin reaction.

P280: Wear protective gloves/protective clothing/eye protection/face protection.

P302+P352: If on skin: Wash with plenty of soap and water.

P333+P313: If skin irritation or rash occurs: Get medical advice/attention.

P362+P364: Take off contaminated clothing and wash it before reuse.

Handling Requirements

- Use of this product should be limited to personnel trained in the techniques of PCR.
- Take normal precautions required for handling all laboratory reagents.

For patient samples only:

- Disinfect spills promptly with 0.5% sodium hypochlorite solution (1:10 v/v bleach) or equivalent disinfectant.

For all reagents:

- Disinfect spills promptly using Microcide SQ. Do not use bleach.
- Regard contaminated materials as biohazardous.
- Wear personal protective apparel, including disposable gloves, throughout the assay procedure. Thoroughly wash hands after removing gloves and dispose of the gloves as biohazardous waste.
- To minimize risk of carryover contamination, regularly change gloves.
- Do NOT:
 - Eat, drink, smoke or apply cosmetics in areas where reagents or samples are handled.
 - Pipette by mouth.
 - Use reagents if reagents appear turbid or cloudy after bringing to specified temperature.
 - Use components beyond expiration date printed on kit label.
 - Interchange vial or bottle caps, as cross-contamination may occur.
- Avoid the use of sharp objects wherever possible.
 - If skin or mucous membrane exposure occurs, immediately wash the area with copious amounts of water; seek medical advice.
- Use all pipetting devices and instruments with care and follow the manufacturers' instructions for calibration and quality control.
- Avoid contamination of reagents and samples.
 - Use aerosol-resistant pipette tips and use a new tip every time a volume is dispensed.
- Do not mix reagents from different kit lots.
- Dispose of hazardous or biologically contaminated materials according to the practices of your institution. Discard all materials in a safe and acceptable manner in compliance with all country, state, local and regulatory requirements.

Storage and Handling

Store the components of this FTD product in its original packaging at -30°C to -10°C. Components are stable until the expiration date as stated on the outer box label.

Freeze product immediately after use. Reagents can sustain up to 10 freeze-thaw cycles.

Specimen Collection and Handling

This section describes the general industry practice for upper respiratory tract specimen handling and storage that will help ensure accurate test results.

This test is for use with extracted nucleic acids from upper respiratory specimens (such as nasal swabs, nasopharyngeal swabs, oropharyngeal swabs, nasopharyngeal wash/aspirate and nasal aspirate) and bronchoalveolar lavage of human origin.

Respiratory pathogen detection depends on the collection of high-quality specimens, their rapid transport to the laboratory, appropriate storage and treatment before laboratory testing. Transport the specimen to the laboratory immediately after collection and process/test as soon as possible. Label all specimens appropriately according to the laboratory's procedure. To protect the viral RNA from degradation, correct specimen handling is very important (as recommended by CDC²).



CAUTION

Handle all samples as if the samples contain potentially infectious agents.

Use universal precautions, such as wearing personal protective apparel, including disposable gloves. Thoroughly wash hands after removing gloves and dispose of the gloves as biohazardous waste.

Prior to sample collection, no special preparation of the patient is required. No pretreatment is required for sample storage.

Collecting the Specimen

Collect all specimens according to the standard technique from the laboratory or clinician. Refer to the CDC Interim Guidelines for Collecting, Handling, and Testing Clinical Specimens from Patients Under Investigation (PUIs) for 2019 Novel Coronavirus.³ Obtain swabs directly from the infection site in order to avoid contamination with surrounding microbiota.

IMPORTANT! Remel M4RT[®] transport medium is not recommended for use with FTD SARS-CoV-2.

Storing and Transporting the Specimen

Nasal, nasopharyngeal and oropharyngeal swabs should be placed immediately in a sterile transport tube containing a viral transport medium or similar method.

Nasopharyngeal wash/aspirate, nasal aspirates (and the non-bacteriostatic saline used to collect these specimens) as well as BAL should be placed immediately into a sterile transport tube or dry container.

Specimens (upper respiratory or BAL) should be delivered promptly to the laboratory and can be refrigerated at temperatures of 2–8°C and/or shipped on ice pack for up to 72 hours after collection, as recommended by CDC.³ Specimens may be frozen to -20°C or ideally -70°C and shipped on dry ice if further delays are expected.

It is important to avoid repeated freezing and thawing of specimens (as recommended by WHO⁴).

NOTES:

- Specimen storage and shipment advice are only recommendations. Check local regulations and institutional policy.
- Pack and label the specimen in compliance with your local and international regulations that cover the transport of clinical samples and etiological agents.

Procedure

Materials Provided

Table 1 details the components for FTD SARS-CoV-2.

Table 1: FTD SARS-CoV-2 Components

Reagent	Composition	Description / Quantity	Storage
SCoV2 PP Mix 96	Synthetic oligonucleotides, buffer	PP mix for SARS-CoV-2 (N gene), SARS-CoV-2 (ORF1ab) and IC 96 reactions: 1 x 144 µL	-30°C to -10°C
SCoV2 PC 96	Double-stranded synthetic DNA molecules, buffer, stabilizing agents	— 96 reactions: 1 x 150 µL	
Negative Ctrl 96	Nuclease-free water	— 96 reactions: 1 x 2000 µL	
Internal Ctrl 96	Double-stranded circular DNA molecules, buffer, <5.0% guanidinium hydrochloride, <0.1% maleic acid	— 96 reactions: 1 x 350 µL	
25x RT-PCR Enz. 96	Enzymes, buffer, glycerol, Ribonuclease (RNase) inhibitors, stabilizing agents	25x RT-PCR Enzyme mix 96 reactions: 1 x 96 µL	
2x RT-PCR Buff. 96	Tris-hydrochloride (Tris-HCl) buffer, Deoxyadenosine triphosphate (dATP), Deoxycytidine triphosphate (dCTP), Deoxyguanosine triphosphate (dGTP), Deoxythymidine triphosphate (dTTP), PCR adjuvants	2x RT-PCR Buffer 96 reactions: 1 x 1200 µL	

Legend: PP = Primer/probe, PC = Positive control, Ctrl = Control, Enz. = Enzyme, Buff. = Buffer

IMPORTANT! The table above reflects the standard kit color scheme. Due to supplier issues during the COVID-19 crisis, individual tube cap colors may be substituted due to availability. Always check the labeling of the reagent prior to use.

Each vial contains additional volume for pipetting inaccuracy. The box and each vial are labeled with a lot number.

Reagents in the kits are sufficient for 96 reactions. Each kit includes IC, NC and PC components.

REF	Contents	Number of Reactions
11416302 (FTD-114-96)	FTD SARS-CoV-2	96

Materials Required but Not Provided

The kit has been validated with the NucliSENS® easyMAG® (bioMérieux)/Applied Biosystems® 7500 Real-Time PCR System (ThermoFisher Scientific), and the VERSANT® kPCR Molecular System (Siemens Healthineers).

The following material and reagents are required for extraction with the NucliSENS® easyMAG®:

Supplier Part Number	Contents
280133	NucliSENS® easyMAG®, Magnetic Silica Beads
280134	NucliSENS® easyMAG®, Lysis Buffer
280130	NucliSENS® easyMAG®, Extraction Buffer 1
280131	NucliSENS® easyMAG®, Extraction Buffer 2
280132	NucliSENS® easyMAG®, Extraction Buffer 3
280135	NucliSENS® easyMAG®, Disposables
N/A	Nuclease-free water

NOTE: Refer to the manufacturer (bioMérieux) for specific part number information.

The following material and reagents are required for extraction with the VERSANT® kPCR Molecular System SP (SP Module):

Supplier Part Number	Contents
04801677	VERSANT® Sample Preparation 1.0 Reagents Box 1
04801685	VERSANT® Sample Preparation 1.0 Reagents Box 2
06691055	96 Deep Well Plate 2 mL
10282930	Disposable Tips 300 µL
10282929	Disposable Tips 1 mL
10489008	Reagent Containers 4 x 200 mL & 2 x 60 mL
N/A	Nuclease-free water

NOTE: Refer to the manufacturer (Siemens Healthineers) for specific part number information.

If viral inactivation using VERSANT® Sample Preparation 1.0 Lysis Buffer is desired prior to extraction with the VERSANT® kPCR Molecular System SP (SP Module), this material is supplied in VERSANT® Sample Preparation 1.0 Reagents Box 1 or available as a separate kit, as follows:

Supplier Part Number	Contents
10634919	VERSANT® Sample Preparation 1.0 Lysis Buffer

NOTE: Refer to the manufacturer (Siemens Healthineers) for specific part number information.

FTD recommends use of an external RNA positive control (RNA PC) with each run, such as:

Supplier Part Number	Contents
0505-0126	AccuPlex™ SARS-CoV-2 Reference Material Kit

NOTE: Refer to the manufacturer (Seracare) for specific part number information.

General Laboratory Equipment and Consumables

- Adjustable micropipette capable of dispensing 1000 µL, 200 µL, 100 µL, 20 µL and 10 µL
- Disposable, aerosol-resistant pipette tips, sterile-packaged
- Disposable, powder-free gloves
- Vortex mixer
- Desktop centrifuge
- Sample rack
- Sample collection devices/material
- Bleach/Microcide (or other product according to laboratory cleaning procedure)
- 96-well PCR plates and plate sealers
- Tubes

Assay Procedures

Workflow Using the NucliSENS® easyMAG® System and the Applied Biosystems® 7500 Real-Time PCR System

Extraction Using the NucliSENS® easyMAG® System

Preparing the Sample

1. Thaw negative control (NC, white cap) and internal control (IC, dark blue cap).
2. Before use, ensure reagents have reached room temperature (15°C to 30°C); mix NC, IC and external RNA PC (by short vortexing) and spin down briefly.
3. Prepare samples for extraction procedure (thaw, if applicable, to room temperature).

Table 2 shows the validated extraction volumes.

Table 2: Validated Extraction Volumes

Type	Volume
Sample volume	200 µL
Elution volume	55 µL

4. Add samples, external RNA PC and NC into the disposables.
5. Program machine accordingly.
6. Select **DISPENSE LYSIS** for lysis buffer to dispense and to start the incubation step. During incubation period, prepare beads as described in the NucliSENS® easyMAG® manual.
7. Once incubation finishes, add 2 µL IC directly to the mix of lysis buffer and sample.
8. Add beads to each well of the disposable and perform extraction protocol.



WARNING

- Never add the IC prior to addition of lysis buffer.
- Never add the IC after extraction.
- Adding IC to each of the samples, the external RNA PC and to the NC is an important step to monitor nucleic acid extraction, as well as the inhibition of nucleic acid amplification.
- Do not extract positive control provided with the kit.

Real-Time PCR Using the Applied Biosystem® 7500

Preparing the Experiment

- Before use, ensure reagents are completely thawed, mixed (by short vortexing) and spun down briefly.

Exception:

- The 25x RT-PCR enzyme must be stored at -30°C to -10°C or on a cooling block at all times.
- Positive Control: Thaw PC and store at room temperature (15°C to 30°C) for 20 to 30 minutes. Vortex PC thoroughly before use.

Table 3: Volume of Reagents Required for 1, 10, 32 and 96 Reactions

Number of Reactions	1	10	32	96
2x RT-PCR Buffer	12.5 µL	125 µL	400 µL	1200 µL
Primer/Probe Mix	1.5 µL	15 µL	48 µL	144 µL
25x RT-PCR Enzyme	1 µL	10 µL	32 µL	96 µL
Total	15 µL	150 µL	480 µL	1440 µL

- Prepare a separate 1.5 mL tube per primer/probe mix and label accordingly. Pipette the required amount of 2x RT-PCR buffer based on the number of reactions (see Table 3).
- Pipette the required amount of SCoV2 PP Mix in the corresponding tube containing 2x RT-PCR buffer (see Table 3).
- Master Mix Preparation:

NOTES:

- In order to obtain accurate volumes and to avoid wasting material, do not immerse the whole tip into the liquid when pipetting the 25x RT-PCR enzyme.
 - Pipette liquid very slowly to prevent air bubbles.
 - Wipe tip against edge of vessel to remove excess liquid outside the tip before dispensing.
 - Change tip after each pipetting step.
- Pipette the required amount of 25x RT-PCR enzyme in each of the tubes containing SCoV2 PP Mix and 2x RT-PCR buffer (see Table 3).
 - Vortex master mix briefly and spin it down.
 - Use master mix immediately and do not store after use.

Prepare a 96-Well Plate for the Applied Biosystems® 7500

NOTE: The PC, NC and external RNA PC must be run on each plate to perform analysis.

NOTE: The RNA PC is recommended but is not provided. Refer to the *Materials Required but Not Provided* section on page 10 for more information.

Refer to Figure 1 for an example of the placement of patient samples and controls.

Figure 1: Samples and Controls – Plate Map Example

	1	2	3	4	5	6	7	8	9	10	11	12
A	Sample 1											
B	Sample 2											
C	Sample 3											
D	Sample 4											
E	Sample 5											
F	RNA PC											
G	PC											
H	NC											

Legend: Green = SCoV2 master mix (A1–H1) • RNA PC = Ribonucleic acid Positive Control (F1)
• PC = Positive Control (G1) • NC = Negative Control (H1)

To prepare a 96-well plate (compatible with the Applied Biosystems® 7500):

1. Pipette 15 µL of the SCoV2 master mix into wells A1 to H1.
2. Add 10 µL of the extracted samples into wells A1 to E1.
3. Add 10 µL of the extracted RNA PC into well F1.
4. Add 10 µL of the PC into well G1.
5. Add 10 µL of the extracted NC into well H1.
6. Seal plate with appropriate adhesive film.
7. Gently vortex plate, then centrifuge briefly.
8. Place plate into the Applied Biosystems® 7500.

NOTE: Refer to manufacturers' operating instructions for use of the Applied Biosystems® 7500.

Program the Thermocycler

Table 4 lists the detection wavelengths for the dyes used in this kit.

Table 4: Detector Programming

SCoV2 PP Mix and Thermocycler Detection Settings		
Pathogen	Dye	Detection Wavelength (nm) ^[a]
SARS-CoV-2	green	520
—	yellow	550
—	orange	610
IC (EAV)	red	670

[a] Detection wavelengths listed are from the Applied Biosystems® 7500. Wavelengths may vary for other thermocyclers.

NOTE: Both targets (N gene and ORF1ab) are labeled with the same dye and are detected in the same channel.

NOTE: Change setting for passive reference dye to **NONE** (by default, ROX dye is selected).

PCR Program

The chart below details the programming steps for the thermocycler.

Stage	Cycles	Acquisition	Temperature	Time
Hold	/	/	50°C	15 minutes
Hold	/	/	94°C	1 minute
Cycling	40	/	94°C	8 seconds
		Yes	60°C	1 minute

For more information on how to program the thermocycler, go to www.fast-trackdiagnostics.com.

Completion of Run

Remove the sealed PCR plate according to the thermocycler manufacturers' instructions. Review the results and discard the PCR plate as biohazardous waste, according to local regulatory requirements.

Workflow Using the VERSANT® kPCR Molecular System

Extraction Using the VERSANT® kPCR Molecular System (SP Module)

Preparing the Sample

1. If viral inactivation is desired prior to extraction, add 350 µL of sample to 350 µL of VERSANT® Sample Preparation 1.0 Lysis Buffer, mix by inversion, then proceed with the automated extraction procedure beginning with Step 2 below.
2. Thaw negative control (NC, white cap) and internal control (IC, dark blue cap).
3. Before use, ensure reagents have reached room temperature (15°C to 30°C); mix NC, IC and external RNA PC (by short vortexing) and spin down briefly.

NOTE: Refer to the external RNA PC manufacturer's preparation instruction. The RNA PC is recommended but is not provided. Refer to the *Materials Required but Not Provided* section on page 10 for more information.

4. Prepare a 1:5 IC dilution in a 2 mL screw cap tube with conical bottom. See Table 5 for required volumes depending on total number of samples per extraction.

Table 5: Required IC Dilution Volume per Extraction

Samples per Extraction*	10	32	64	81	96**
Internal Control	48 µL	96 µL	166 µL	201 µL	234 µL
Nuclease-free water	192 µL	384 µL	664 µL	804 µL	936 µL
Total IC Dilution Volume	240 µL	480 µL	830 µL	1005 µL	1170 µL

* including external RNA positive control (RNA PC) and negative control (NC).

** When using Dynamic Assay Preparation FTD1 (SP1.0 500-15-10), one FTD-SARS-CoV-2 kit is enough to run 82 PCR reactions (81 extracted samples [including RNA PC and NC] plus one non-extracted positive control).

The IC input volume for the Internal Control Mix varies based on the number of samples being run on the SP Module. Calculate the IC input volume using the following formula:

Formula	(# of Samples x IC output per sample + [5 x PreAsp volume] x 2) x (1 + safety margin) + 15 µL	= IC input volume
Example for 15 samples	(15 x 10 µL + [5 x 10 µL] x 2) x (1 + 10%) + 15 µL = 290 µL	

Legend: PreAsp volume = Output volume to Deep well plate • Safety margin = 10% • 15 µL = Dead volume

5. Prepare samples for extraction procedure (thaw, if applicable, to room temperature).
6. If necessary, transfer 650 µL of each sample into sterile VERSANT® kPCR SP Module compatible polypropylene tubes that meet the size requirements for the sample carrier.

NOTE: Sample tubes must meet the following requirements:


- Do not use sample tubes with internal threads.
- For the 32-position sample carrier, only use 11 x 60 mm to 14 x 100 mm sample tubes with an inner diameter greater than or equal to 9 mm.

7. Place barcode labels on all sample tubes. The system requires that all samples are uniquely identified with a barcode. See Figure 2 for guidance on placing the barcode on the sample tube.

NOTES:

- For a given run, ensure that each sample barcode is unique.
- Position the barcode within a range of 20 to 100 mm (0.79 to 3.94 in) from the bottom of the sample tube.
- Affix the barcode tightly to the sample tube in a vertical position.

Figure 2: Placement of the Barcode Label on the Sample Tube

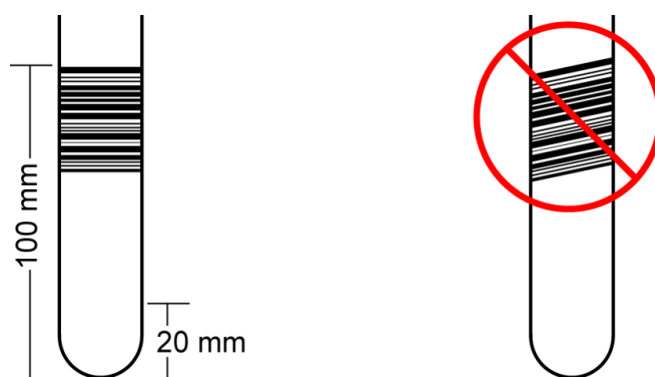

8. Load the samples onto the sample carriers.

NOTE: Ensure that the barcode on each of the sample tubes is positioned correctly so that it is clearly visible through the opening in the carrier.

9. When placing the samples in the sample carrier, carefully remove the tube caps.

- To reduce risk of cross-contamination, discard sample tube caps as you remove them.
- Use new caps when recapping the tubes.

NOTE: If any splashing of samples occurs, remove and discard your gloves and put on clean ones.

10. To load the Sample Carriers onto the autoload tray, insert each carrier into the tracks along the front and rear carrier guides. Push the carrier in until it touches the stop clips just beyond the autoload tray. Do not attempt to force the carrier past the stop clips.

NOTES:

- Up to four sample carriers can be loaded in positions 20–23 on the autoload tray.
- Place the Sample Carriers in positions 20–23 on the autoload tray.
- Ensure that the barcodes on the Sample Carriers are facing to the right.
- Before starting the run, ensure that all the carriers are loaded in the correct positions on the autoload tray.


11. After loading the prepared carriers onto the autoload tray, start the run.

Entering Lot Information

To enter lot information, perform the following actions:

1. Select the **Protocol** tab at the main screen (see Figure 3).


Figure 3: Protocol Tab

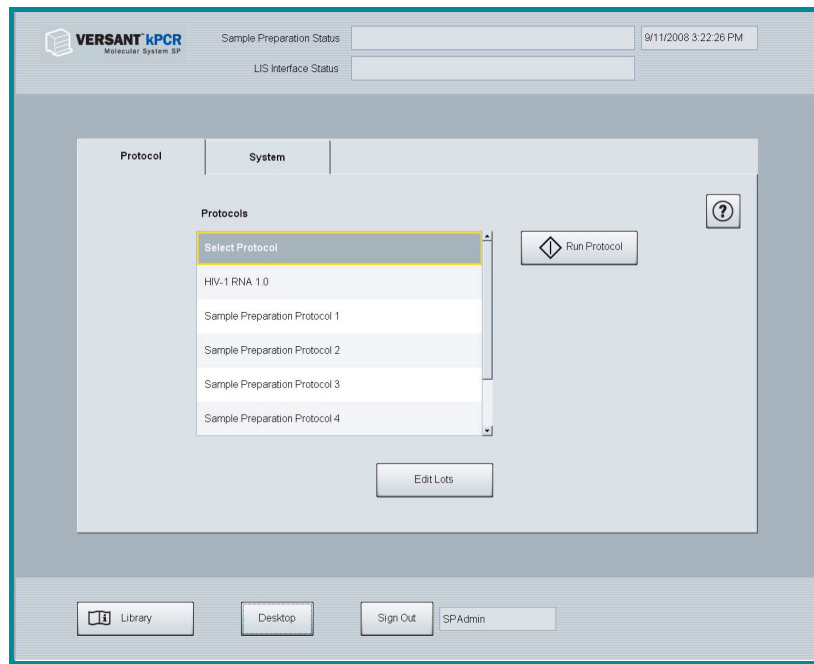

2. If you are using a kit from a new lot of sample preparation reagents, enter the lot-specific information using the hand-held scanner:
 - a. At the **Protocol** tab, select **Edit Lots**.

The **Kit Lot List** window displays (see Figure 4).

Figure 4: Kit Lot Window

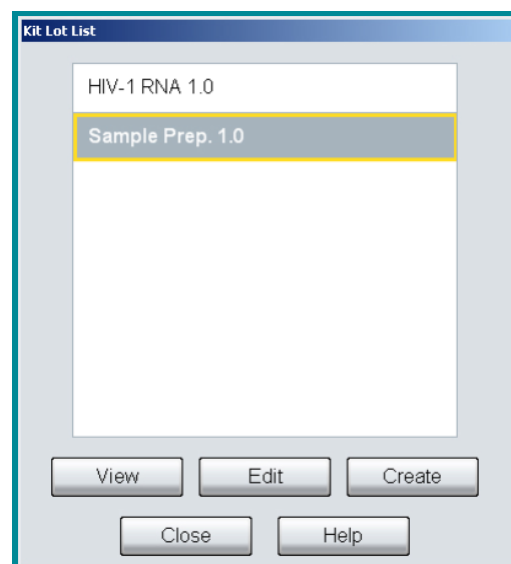

- b. Select **Sample Prep 1.0** from the list.
- c. Select **Create**.

The **Enter Lot Information** window displays (see Figure 5).

Figure 5: Enter Lot Information Window

Assay Code	
Build Number	
Kit Lot	
Expiration (yyyy-mm-dd)	

Save Cancel

- d. Place the cursor in the first row of the right column.
 - e. Scan the 2D barcode that is printed on the assay product insert supplement.
3. When you scan the 2D barcode, the form fields populate automatically. Verify the lot information entered:

NOTE: Scanning the 2D barcode is only required for new lots of reagent kits. However, as a good practice, check the assay kit information before each run:

- a. At the protocol tab, select **Edit Lots**.
 - b. Select **Sample Prep 1.0** from the list.
 - c. Select **View**.
 - d. Verify that the information is correct for the assay kit lot you plan to run.
 - e. Select **OK**.
4. The system returns to the **Kit Lot List** window.
5. If the lot information is correct, select **Close** to return to the main system window and start your run.


NOTE: If assay kit information is not correct, select **Create** and scan the 2D barcode as described above.

Starting the Run

To start the run, perform the following actions:

1. Select the **Protocol** tab at the main screen (see Figure 6).

Figure 6: Protocol Tab

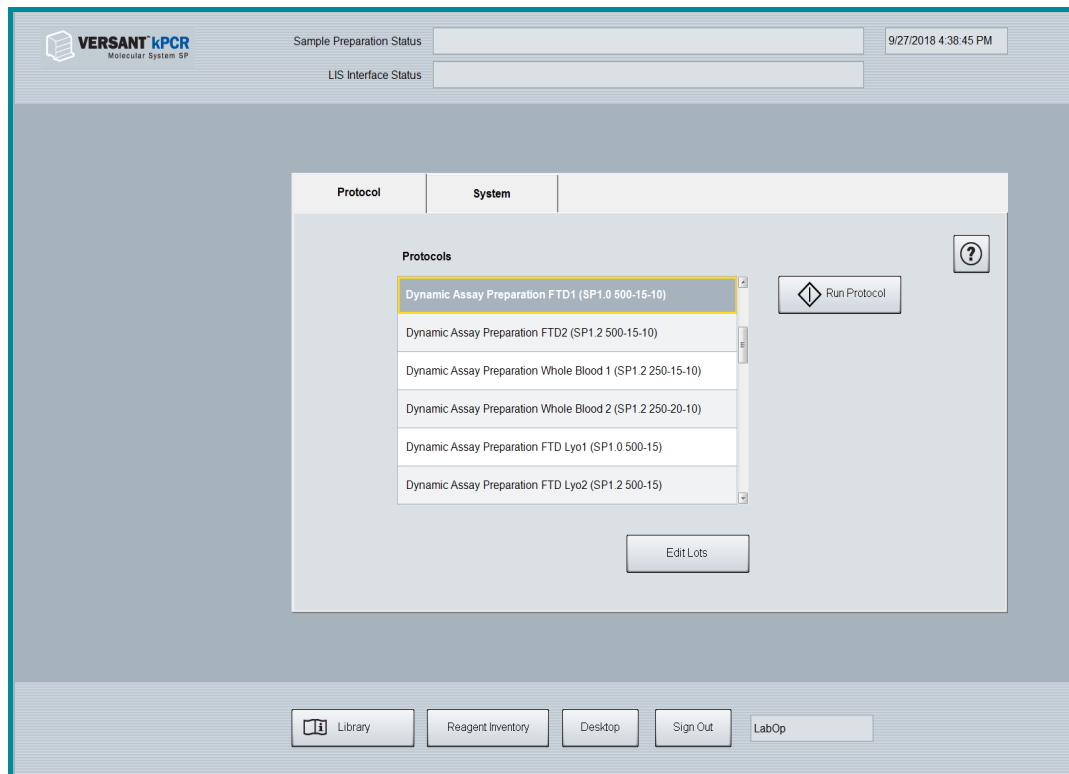

2. Select **Dynamic Assay Preparation FTD1 (SP1.0 500-15-10)**.
3. Click **Run Protocol**.

The SP Module verifies that it can initialize the iSWAP, the waste sensors, and the heater/shaker. If the system is unable to complete initialization, an error window displays.

Refer to the *VERSANT® kPCR Molecular System Sample Preparation (SP) Application Guide⁵* for more information about these errors.

4. After the system completes initialization, you are prompted to empty the solid waste container.

Refer to the *VERSANT® kPCR Molecular System Sample Preparation (SP) Application Guide⁵* for instructions on emptying and reinstalling the solid waste container.

5. When the solid waste container is reinstalled correctly, the **Install the Solid Waste Container** window closes, and the **Dynamic Assay Run Setup** window displays.

Configuring a Dynamic Assay Run

When starting a run, the **Dynamic Assay Run Setup** window opens automatically and displays the information that was saved from the last run (see Figure 7). Use this window to configure a Dynamic Assay run.

Figure 7: Dynamic Assay Run Setup

1. In the **Number of Assays** list, select **1**.
2. Select the **Dispense IC** checkbox.

If you select the checkbox, Internal Control (IC) (supplied as Internal Ctrl 96 in the FTD SARS-CoV-2 kit) is dispensed to all samples and controls extracted by the SP Module.



WARNING

- Never manually add IC to the sample, before or after extraction.
- Adding IC to each of the samples and to the NC is an important step to monitor nucleic acid extraction, as well as the inhibition of nucleic acid amplification.

3. Select the **Universal Negative Control** checkbox.

If you select the checkbox, **Universal Negative Control**, a full process negative control (supplied as Negative Ctrl 96 in the FTD SARS-CoV-2 kit) is extracted along with the clinical samples on the SP Module.

-
4. In the **Amplification Detection Protocol Name** list, type or select **FTD-SARS-CoV-2**.
-

CAUTION

The Amplification Detection Protocol Name on the SP module must match the kPCR protocol name saved on the AD Module to link the Dynamic Assay Protocol information to the AD.

5. Under **Assay 1**, type or select **FTD-SARS-CoV-2**.

If you select a stored Assay name, the **Controls**, **Kit Lot Number**, and **Expiration Date** fields populate with the stored assay information, and the **Dispense IC** and **Universal Negative Control** fields display values from the last saved run.

CAUTION

If you are using the LIS, the assay name must match the dye channel name entered in the AD Module set up.

6. Under **Number of Calibrator/Control(s)**, select:
 - a. **1** under **Control (No SP)**. This PCR-Pos control (supplied as SCoV2 PC 96 in the FTD SARS-CoV-2 kit) will be added directly to the PCR plate after the sample extraction process is complete.
 - b. **1** under **Control (SP)** (External RNA positive Control [RNA PC]).
7. (Optional) Under **Kit Lot Number**, type the kit lot number.
8. Under **Expiration Date**, select the expiration date for the lot.

Expired reagents are highlighted in red on the **Dynamic Assay Run Setup** window and the VERSANT kPCR Sample Preparation Run Report. The system does not enforce the expiration of the kit lot.

9. Under **Expiration Date**, select the expiration date for the lot.
10. Click **Continue**.

Loading Samples

1. Ensure that the barcode on each of the samples is visible through the opening of the sample carrier. The barcode of each of the samples is read as the carriers are loaded.
2. Ensure that you have loaded all sample carriers on the autoloader tray.
3. Click **Load Samples**.

The system begins loading the sample carriers from the autoloader tray onto the deck and reads all carrier IDs, sample barcodes and sample slot locations for sample tubes contained on the sample carrier.

If errors occur, refer to the *VERSANT® kPCR Molecular System Sample Preparation (SP) Application Guide*⁵ for more information about manually entering barcodes or troubleshooting the error.

When all sample carriers are loaded, the **Loading Samples Confirmation** window displays, indicating the number of samples and carriers loaded.

- If the information in the **Loading Sample Confirmation** window correctly indicates the number of samples and sample carriers you loaded, click **Yes**.

The **Work List Editor** displays and shows the PCR Plate Map and the Work List (see Figure 8). The Work List includes the sample IDs of all samples loaded on the SP instrument. If the system is connected to the LIS, all test orders will be displayed in the Work List.

Figure 8: Work List Editor

The screenshot shows the 'Work List Editor' window. At the top, there is a table titled 'FTD-SARS-CoV-2' with columns 1 through 12. The rows are labeled A through H. The first column (A) contains sample IDs: A1 (FTD-SAR...), A2 (Universal...), A3 (FTD-SAR...), A4 (100001), A5 (100002), A6 (100003), A7 (100004), A8 (100005), A9 (100006), A10 (100007), A11 (100008), A12 (100009). The other columns (B-H) contain sample IDs for each row. Below the table, there is a section titled 'Total Samples: 78' with buttons for 'No SP Analyte Reagent', 'SP Analyte Reagent', and 'Sample ID'. Below this is a table with 'Sample ID' and 'FTD-SARS-CoV-2' columns. The 'Sample ID' column lists 10 sample IDs (100001 to 100010). The 'FTD-SARS-CoV-2' column has checkboxes, all of which are checked. At the bottom, there are buttons for 'Load/Refresh', 'Clear All/Selected', 'Check All/Selected', 'Cancel', and 'Continue'. There is also a status bar with 'Ordered', 'LIS Requested, Not Ordered', and 'Not Ordered'.

Legend: A1 – Positive Control • B1 – Extracted Negative Control (NC) • C1 – Extracted external RNA Positive Control (RNA PC)


- Select all samples in the work list editor (right column) and click **Check All/Selected**.
- Click **Continue**.

About the PCR Plate Map

The PCR Plate Map is populated with the control and sample IDs based on the Dynamic Assay Protocol Setup and the Work List selections in the LIS.

The PCR plate map is read-only, except for the names of the controls. You can modify the control well names to match the reagents that you are loading. The control name must be between 2 and 20 letters, numbers, and symbols and is case sensitive. The control name becomes the sample ID. The Universal Negative Control name cannot be changed.

The PCR plate map reflects the order that the controls and samples are loaded into the PCR plate. This is based on the number of SP and No SP controls selected for each assay. The plate is loaded in the following order:

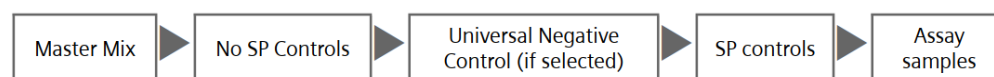

About the Work List

The first column of the work list contains the sample IDs for all samples loaded on the system starting with position 1 of the sample rack loaded in track position 23 on the deck. The adjacent column headers include the assay names that were defined in the **Dynamic Assay Run Setup** window.

Test order information from the LIS:

- If a sample ID is associated with an assay in the test order from the LIS, then the assay cell is checked, otherwise it is unchecked.
- If a sample ID from the LIS cannot be accommodated on the PCR Plate due to a lack of available wells for an assay, each associated check box for that sample will be disabled and highlighted red.
- A sample ID with no assays selected will be in highlighted in red.

The following table describes the actions you can perform for any sample in the work list (including samples with test orders from the LIS).


To	Do
Select the required assay(s) for each sample.	Select the checkbox in the assay column next to the sample ID.
Reload the work list with the initial sample information from the LIS.	Click Load/Refresh .
Clear all checkboxes in the worklist.	Click Clear All/Selected .
Select all checkboxes in the worklist.	Click Check All/Selected .
Abort the workflow, unload the racks from the deck to the autoloader tray, and return to the Dynamic Assay Run Set Up window.	Click Cancel .
Continue the workflow from the SP Control Carrier Loading window and load the consumables.	Click Continue .

Note that the PCR plate map entries display as you select items in the work list editor.

Loading Consumables and Reagents

After loading samples, click **Continue** in the **Worklist Editor** window, the **Tip Trays** window displays (see Figure 9).

Figure 9: Work List Editor

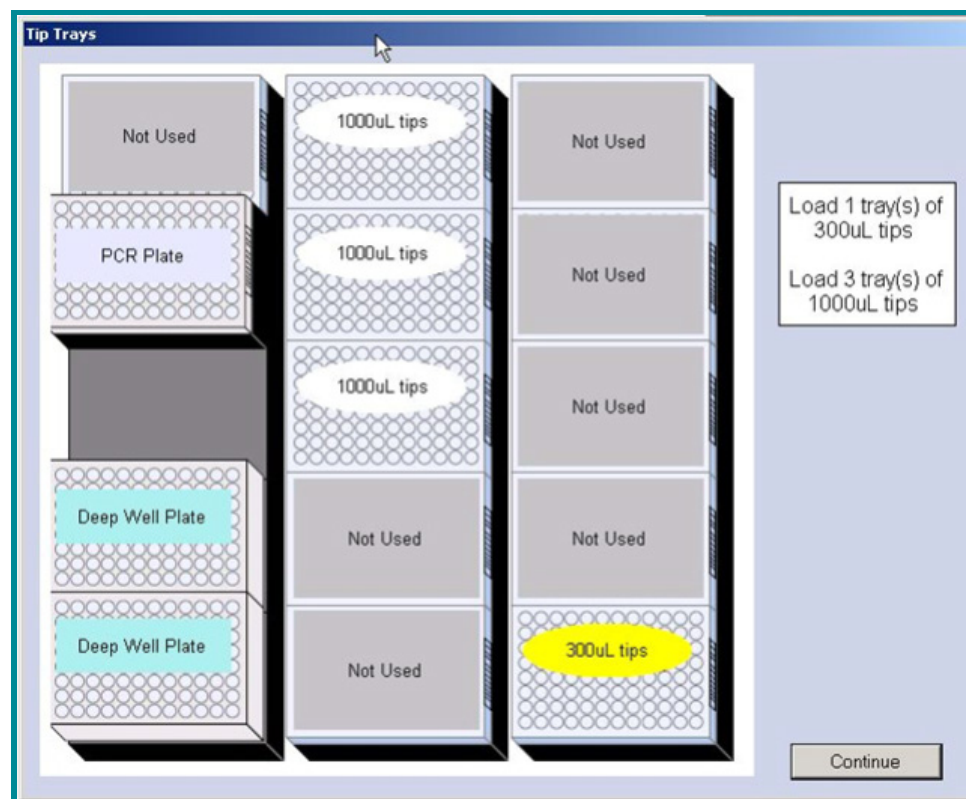

1. Fill the Carriers with tips and plates.

The **Tip Tray** window displays the locations for 300 μ L and 1000 μ L tips-racks deep well plates and PCR plate based on the number of assays, controls and samples selected.

2. Load the Tip Tray Tip Carriers onto the deck.
3. Prepare the VERSANT SP 1.0 Reagents.
 - a. Visually inspect each reagent for the presence of solids. If small amounts of solids are present, allow them to settle before pouring the reagent. Ensure the solids do not go into the troughs.
 - b. Mix bottles containing Lysis Buffer, Wash Buffer 1, Wash Buffer 2, and Wash Buffer 3 by gently inverting 10 times.

NOTE: Slight color changes may occur with the Lysis Buffer. These changes do not indicate a change in the quality of the buffer.

- c. Pour the entire volume of each reagent into a new, properly labeled, 200 mL trough. Minimize formation of bubbles by pouring along the inside wall of the trough.

4. Load the Large Trough Carrier:

- Lysis Buffer in Position 1 (nearest the carrier barcode)
- Wash 1 in Position 2
- Wash 2 in Position 3
- Wash 3 in Position 4 (nearest the carrier handle)



5. Load the Small Trough and Control Carriers:

- Mix the Magnetic Beads by vortexing the bottle for 60 seconds.
- Pour Magnetic Beads into a new, properly labeled, 60 mL trough.
- Mix the Elution Buffer by gently inverting bottle 10 times.
- Pour the Elution Buffer into a new, properly labeled, 60 mL trough.
- Place Magnetic Beads in Position 4 and Elution Buffer in Position 5 in the small trough carrier as follows:



- Place the small trough carrier at deck Position 18 on the autoload tray of the SP Module.

6. Click **Continue**.

The **SP Analyte Reagent Carrier Loading** window displays (see Figure 10).

Figure 10: SP Analyte Reagent Carrier Loading Window

SP Analyte Reagent Carrier Loading

Ensure caps are removed from all the tubes and select Load Consumables

NOW: Prepare and load the SP analyte reagent carrier

SP Reagents	
D-SARS-CoV-2	1 FTD-SARS-CoV-2-C1
	2
	3
	4
	5 DO NOT USE
	6
	7
	8
	9
	10 DO NOT USE
	11
	12
	13
	14
	15 DO NOT USE
	16
	17
	18
	19
	20 DO NOT USE
	21
	22
	23
	24
	25 DO NOT USE
	26
	27
	28
	29
	30 Universal Negative Control
	31 Internal Control (191 µL)
	32 PK

Duration of display (hrs:mins:secs):

00:00:10

LATER: Prepare the No SP analyte reagent carrier for loading after the SP steps complete.

No SP Reagents	
SARS-CoV-2	1 MMX (169 µL)
	2 FTD-SARS-CoV-2-A1
	3
	4
	5
	6
	7
	8
	9
	10
	11
	12
	13
	14
	15
	16
	17
	18
	19
	20
	21
	22
	23
	24
	25
	26
	27
	28
	29
	30
	31
	32

7. Add SP reagents as indicated:

- a. Take out the barcode label (0007) from the SP1.0 Reagent Box 2 (Proteinase K bag) and place it on a new, clean 5 mL tube (provided in the same bag).
- b. Mix Proteinase K by gently inverting the bottle 10 times.
- c. Pour the entire volume of Proteinase K into the barcoded 5 mL tube.
- d. Place the barcoded 5 mL tube into Position 32 of the control carrier. Ensure that the barcode is positioned vertically so that it is clearly visible through the opening in the carrier.
- e. Place the internal control dilution into Position 31.
- f. Place the Universe Negative Control (provided as Negative Ctrl 96) into Position 30.
- g. Place the external RNA PC into Position 1 (FTD-SARS-CoV-2-C1).

8. Load control carrier onto deck Position 19 on the autoloader tray of SP Module.

9. Click **Print** to print a copy of the SP/No SP Analyte Reagent Carrier Map Report.

10. Ensure that caps are removed from all tubes and troughs and click **Load Consumables**.

The system performs the following activities:


- Checks that the consumable carriers (Tip Carrier, Plate Carrier, Control Carrier, and Sample Preparation carriers) are available on the autoloader tray.
- Loads the Tip Carrier.
- Loads the Sample Preparation Reagents Carrier and confirms adequate volume is present.
- Loads the Control Carrier.
- Confirms adequate volume is present in the Controls and IC.
- Confirms adequate volume is present for PK.

About Liquid Volume Levels

If the system does not confirm adequate volume is available in the PK, IC or master mix, or if it does not detect liquid in the Controls (No SP):

- The system ejects the Control Carrier from the deck.
- The **SP Analyte Reagent Carrier Loading** window displays the **Liquid Level Too Low** warning (see Figure 11).
- The corresponding reagent position is highlighted.

Figure 11: SP Analyte Reagent Carrier Loading Window with Liquid Level Warning

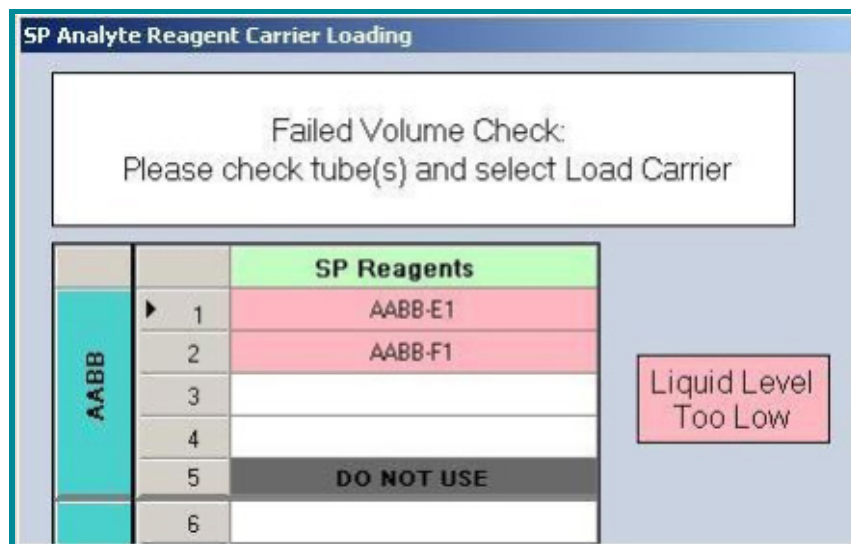


To continue the run, fill the control with adequate volume and then click **Load Consumables**.

After all consumables load, a window displays a timer that indicates when the Control (No SP) Reagents will need to be loaded.

The system begins the sample extraction process.

When all samples have been transferred to the deep well plate, the **Sample Tube Management** window displays the following message (see Figure 12): Sample transfer is complete. Ordered sample tubes are highlighted in green.

Figure 12: Sample Tube Management Window

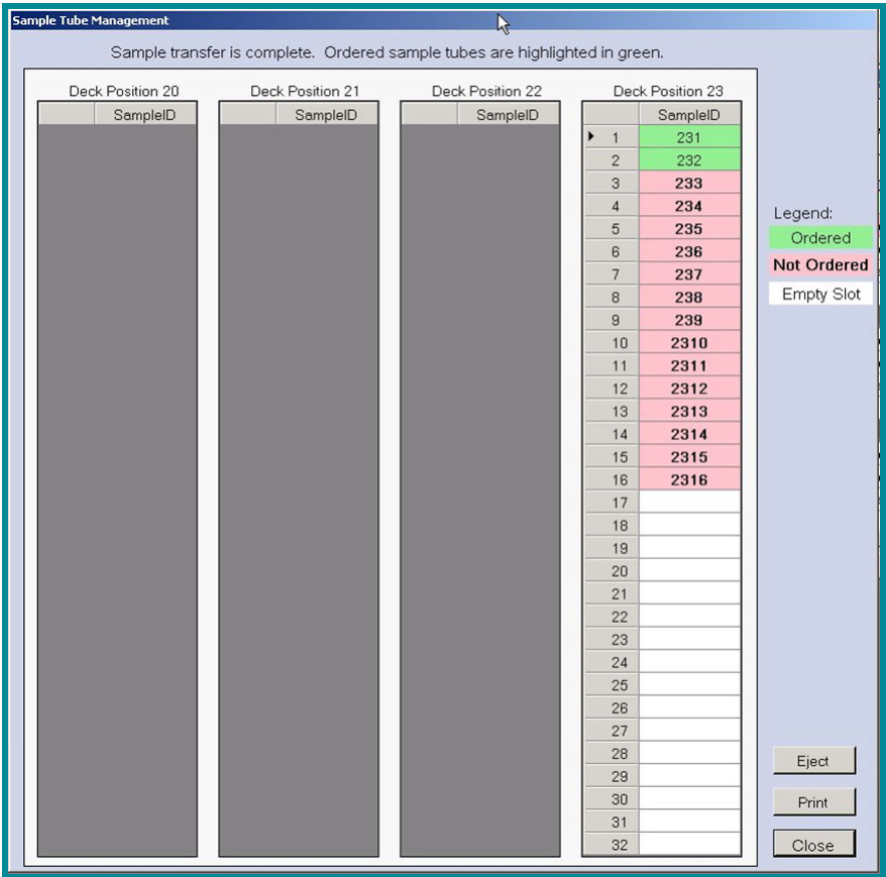

The window remains open for 15 minutes. If you do not select one of the following options, the window automatically closes and the carriers are not ejected.

Option	Description
Print	Print the Sample Tube Management window.
Close	Close the window.
Eject	Eject the sample carriers from the deck.

Preparation of Master Mix and Positive Control (No SP Control)

NOTE: Prepare the master mix and the Control Reagents (No SP) right before the end of the extraction.

- Before use, ensure reagents are completely thawed, mixed (by short vortexing) and spun down briefly.

Exception:

- The 25x RT-PCR enzyme must be stored at -30°C to -10°C or on a cooling block at all times.
- Positive Control: Thaw PC and store at room temperature (15°C to 30°C) for 20 to 30 minutes. Vortex PC thoroughly before use.

- Before use, ensure reagents are completely thawed, mixed (by short vortexing) and spun down briefly.
- Prepare a separate 2 mL tube for the SCoV2 master mix and label accordingly. Pipette the required amount of 2x RT-PCR buffer, SCoV2 Primer/Probe mix and 25x RT PCR Enzyme based on the number of reaction (see Table 6).

Table 6: Volume of Reagents Required for 32, 64, 82 and 96 Reactions

Number of Reactions	32	64	82	96*
2x RT-PCR Buffer	507.5 µL	947.5 µL	1195 µL	1387.5 µL
Primer/Probe Mix	60.9 µL	113.7 µL	143.4 µL	166.5 µL
25x RT-PCR Enzyme	40.6 µL	75.8 µL	95.6 µL	111 µL
Total	609 µL	1137 µL	1434 µL	1665 µL

* When using Dynamic Assay Preparation FTD1 (SP1.0 500-15-10), one FTD-SARS-CoV-2 kit is enough to run 82 PCR reactions (including one No SP control).

- The input volume for the master mix varies based on the number of samples being run on the SP Module. Calculate the input volumes using the following formula:

Formula	$(\# \text{ of samples} \times \text{master mix per well} + [2 \times \text{PreAsp volume}] \times 2) \times (1 + \text{safety margin}) + 15 \mu\text{L}$	= master mix input volume
Example for 48 samples	$(48 \times 15 \mu\text{L} + [2 \times 15 \mu\text{L}] \times 2) \times (1 + 10\%) + 15 \mu\text{L} = 873 \mu\text{L}$	

Legend: PreAsp volume = Output volume to PCR plate • Safety margin = 10% • 15 µL = Dead volume

- After calculating the master mix input volume, calculate the volumes of the respective master mix components:

- 2x RT-PCR Buffer volume = $12.5/15 \times 873 \mu\text{L} = 727.5 \mu\text{L}$
- Primer/Probe Mix volume = $1.5/15 \times 873 \mu\text{L} = 87.3 \mu\text{L}$
- 25x RT-PCR Enzyme = $1/15 \times 873 \mu\text{L} = 58.2 \mu\text{L}$

NOTES:

- To obtain accurate volumes and avoid wasting material, do not immerse the whole tip into the liquid when pipetting the 25x RT-PCR enzyme.
 - Pipette liquid very slowly to prevent air bubbles.
 - Wipe tip against edge of vessel to remove excess liquid outside the tip before dispensing.
 - Change tip after each pipetting step.
4. Vortex master mix briefly and spin it down.
 5. Use master mix immediately and do not store after use.

Loading the Control (No SP) Carrier

**WARNING**

When you see the system prompt, immediately replace the Controls (SP) reagents in the Control carrier with the Controls (No SP) and the master mix reagents.

1. Click **Eject Carrier** to eject the Control carrier from the instrument.
2. Remove the Control carrier from the autoloader tray, remove the reagents from the carrier, and replace them with the master mix and Controls (No SP) reagents.
3. Click **Load Consumables**.

The system performs the following actions:

- Loads the Control Carrier.
- Confirms adequate liquid is present in both the master mix and in the Controls.

Completing the SP Module Run

When a sample preparation run completes, the SP Software audibly beeps and displays the **Run Complete** window.

- If the **Run Complete** dialog box indicates that sample preparation errors occurred during the run, note the error message and complete the run. Refer to the run report to determine where the sample preparation errors occurred.
- If the **Run Complete** dialog box indicates that a plate map file could not be generated, you cannot complete the run. Refer to the *SP Module Troubleshooting* section in the *VERSANT® kPCR Molecular System Sample Preparation (SP) Application Guide*⁵ for more information about possible errors and solutions.
- When the **Run Complete** dialog box displays without errors, complete the run.

To complete the run:

1. Verify that the autoloader tray is clear and perform one of the following actions:
2. After the run completes, immediately remove the PCR plates from the SP Module, immediately seal the PCR plate with an AD Module compatible cover, vortex on high for 10 seconds, plate-centrifuge 1,000xg/2 minutes and start the run on the AD Module.
3. Clean the SP Module.

Real-Time PCR Using the VERSANT® kPCR Module System AD (AD Module)

The VERSANT® kPCR Molecular System AD is a QuantStudio™ 5 Dx real-time thermocycler (originally supplied by ThermoFisher) modified to run the MiPLX Software Solution (Siemens Healthineers).

Running a PCR on the AD Module

To run a kPCR experiment based on a Dynamic Assay protocol, transport the sealed PCR plate to the AD Module.

Refer to the *Steps for Loading and Unloading the Sealed PCR Plate* section in the *VERSANT® kPCR AD Application Guide (ThermoFisher Scientific QuantStudio 5 Dx)*.⁶

Importing the Plate Map

During sample preparation on the SP Module, the system generates a plate map that indicates the position of all the samples along with their corresponding sample IDs. When you start a run on the AD Module, you must first transfer the plate map from the SP Module (the plate map may be found on the SP controlling computer C:\VersantkPCR\Platemap. Example: H990KJYC.map) and to a folder on the PCR computer running the MiPLX software (e.g., C:\MapFiles).

To import the plate map:

1. Open the MiPLX software.
2. Select **File** and then **Import Plate Map**.
3. Enter the run id. The run id is the barcode number on the PCR plate (example: H990KJYC).
 - The plate map for the PCR plate imports into the VERSANT® kPCR Molecular System AD Software.
 - The parameters for the run are automatically selected from the SP Module based on the selected SP Amplification Detection Protocol Name: FTD-SARS-CoV-2.
4. Click **Add** to continue the run or scan the barcode on the PCR plate with a barcode reader to import the plate map.

The plate map automatically displays. Based on the plate map created on the SP Module, the system automatically sets the Sample IDs for all wells.

The Sample IDs for all wells appear on the **Full Screen Plate** button. You can view the complete Sample ID for a well by clicking the individual well and then clicking **Well Information**.



WARNING

Do not change the operating system settings that were configured by the Siemens service representative. This can change the character set available for Sample IDs and result in the use of invalid characters.

Do not edit the Sample IDs associated with the imported plate map file. The Sample ID used for the SP Module run must remain the same for the AD Module run. The Sample ID on Printed Reports may be truncated and LIS errors may occur if the Sample ID is revised.

- The system identifies each well as Unknown.
 - If applicable, the system indicates any errors that occurred during sample preparation on the SP Module.
-

Entering Assay Kit Information

For Dynamic Assay Protocols, after importing a plate map into the AD Module, use the Assay Kit Information feature to enter assay reagent information. Otherwise, the system assumes that the assay reagent kit lot information is the same as the sample preparation reagent kit lot information (such as Kit Lot Number and Expiration Date).

Starting a Run on the AD Module



WARNING

Before starting a new run on the AD Module, select **File** and then **Save As**, and enter a name for the run (type a new name).

The name cannot match the name of any existing Dynamic Assay Protocols. Saving the run with the name of an existing protocol will corrupt that Amplification/Detection protocol template.

Refer to the *Starting the Run* section in the *VERSANT® kPCR AD Application Guide (ThermoFisher Scientific QuantStudio 5 Dx)*.⁶

Completing the AD Module

You can review the results versus the Sample Preparation Run Report to detect any results that were not cleared/not selected from data gathering due to sample preparation errors.

To review the results:

1. Click **Analysis** on the toolbar.
The **Analysis** pane displays.
2. In the **Analysis** pane, click **Results**.
3. Under **Area to Analyze**, click **Text Report** to view the results.
4. After you review the results, select **File** and then **Close** to close the protocol.



WARNING

Analysis term settings that are entered manually on the AD module are not saved. The user must manually enter selected analysis term settings each time an analysis result is generated. Analysis term settings revert to the manufacturer's default settings after each analysis is generated. To prevent reporting incorrect patient results, confirm that the appropriate analysis term settings were used.

Exporting Data to the Laboratory Information System (LIS)

If your system is connected to a LIS, you can export your kPCR protocol data directly to the LIS.

Refer to the *Exporting Data to (LIS) Laboratory Information System* section in the *VERSANT® kPCR AD Application Guide (ThermoFisher Scientific QuantStudio 5 Dx)*.⁶

Quality Control

The FTD SARS-CoV-2 test includes the following controls:

- Negative Control (NC)
- Positive Control (PC)
- Internal Control (IC)

In addition, FTD recommends the use of an external RNA positive control (RNA PC) such as the AccuPlex™ SARS-CoV-2 Reference Material Kit (Material number 0505-0126).

The assay uses equine arteritis virus (EAV) as an IC, which is introduced into each sample and the negative control (NC) during the extraction process. The IC is extracted, processed and amplified simultaneously with each sample in order to monitor the extraction process and to allow the identification of PCR inhibition.

The NC is processed as a sample (extraction and RT-PCR). It confirms the absence of contamination.

The FTD SARS-CoV-2 kit contains a positive control (PC), which is added to each RT-PCR run. It monitors the RT-PCR process and performance of the primers and probes.

The external RNA PC is processed as a sample (extraction and RT-PCR). It confirms the validity of the extraction, reverse transcription and amplification. To obtain a low positive control, the content of the positive vial should be diluted to obtain a final concentration of 700 copies per milliliter (cop/mL).

Be aware that the AccuPlex™ SARS-CoV-2 Reference Material must be extracted prior to RT-PCR according to the instructions from the *Extraction Using the NucliSENS® easyMAG® System* section on page 11 and the *Extraction Using the VERSANT® kPCR Molecular System* section on page 15. The AccuPlex™ SARS-CoV-2 Reference Material Kit content must not substitute the positive control (PC) and negative control (NC) reagents provided with the SARS-CoV-2 kit.

Criteria for a Valid Run

The run is considered valid and patient results are reported if all the following conditions are met:

1. NC shall not show any amplification traces other than the one for the IC. The IC must fall below a Ct of 33. Manually inspect the NC for unspecific amplification detected in the green detection channel. If there is a potential contamination (appearance of a curve in the green detection channel), results obtained are not interpretable and the whole run (including extraction) must be repeated.
2. PC must show a positive (*i.e.*, exponential) amplification trace for SARS-CoV-2. The PC must fall below a Ct of 33.
3. All samples, that do not show SARS-CoV-2 amplification, must show a positive amplification trace for the IC with a Ct less than 33. If IC is negative or shows a Ct value greater than 33 in a SARS-CoV-2 negative sample, the result is invalid. If IC is negative or shows a Ct value greater than 33 in a SARS-CoV-2 positive sample, the run is invalid.
4. External RNA PC must show a positive (*i.e.*, exponential) amplification trace for SARS-CoV-2.

Results

Interpretation of Results

Table 7 details the possible results with FTD SARS-CoV-2.

Table 7: Result Interpretation of Clinical Samples and Controls

Sample/ Control	SARS-CoV-2	IC	Overall Result	Interpretation
Patient Sample	Negative	Ct < 33	Valid	SARS-CoV-2 not detected.
	Negative	Ct ≥ 33	Invalid	There was an error during extraction/PCR. Retest the sample.
	Negative	Not detected	Invalid	There was an error during extraction/PCR. Retest the sample.
	Positive	Ct < 33	Valid	SARS-CoV-2 detected.
	Positive	Ct > 33	Valid	SARS-CoV-2 detected.
	Positive	Not detected	Valid	SARS-CoV-2 detected.
NC	Negative	Ct < 33	Valid	Run is valid.
	Negative	Ct ≥ 33	Invalid	There was an error during extraction/PCR. Run is invalid.
	Negative	Not detected	Invalid	There was an error during extraction/PCR. Run is invalid.
PC	Ct < 33	Not applicable	Valid	Run is valid.
	Ct ≥ 33	Not applicable	Invalid	There was an error during PCR. Run is invalid.
	Not detected	Not applicable	Invalid	There was an error during PCR. Run is invalid.

The results will be reported as a cycle threshold (Ct) unit.


If criteria listed in the *Criteria for a Valid Run* section are met, any patient sample displaying an exponential trace shall be considered as positive for the pathogen targeted by the kit. An absence of an exponential trace indicates an absence or undetectable load of nucleic acid.

The IC must be positive in each extracted material that is not positive for SARS-CoV-2.

IMPORTANT! Pay attention to the section below for important information regarding baseline settings and multicomponent plots.

Baseline Setting and Multicomponent Plot

The amplification curve baseline is one of the parameters that can affect PCR results. In case the baseline is incorrectly set, a Ct value can be displayed even if no real amplification occurred. Figure 13 illustrates the difference between a real amplification (A) and an incorrect baseline setting conveyed with a Ct value even if no amplification occurred (B).

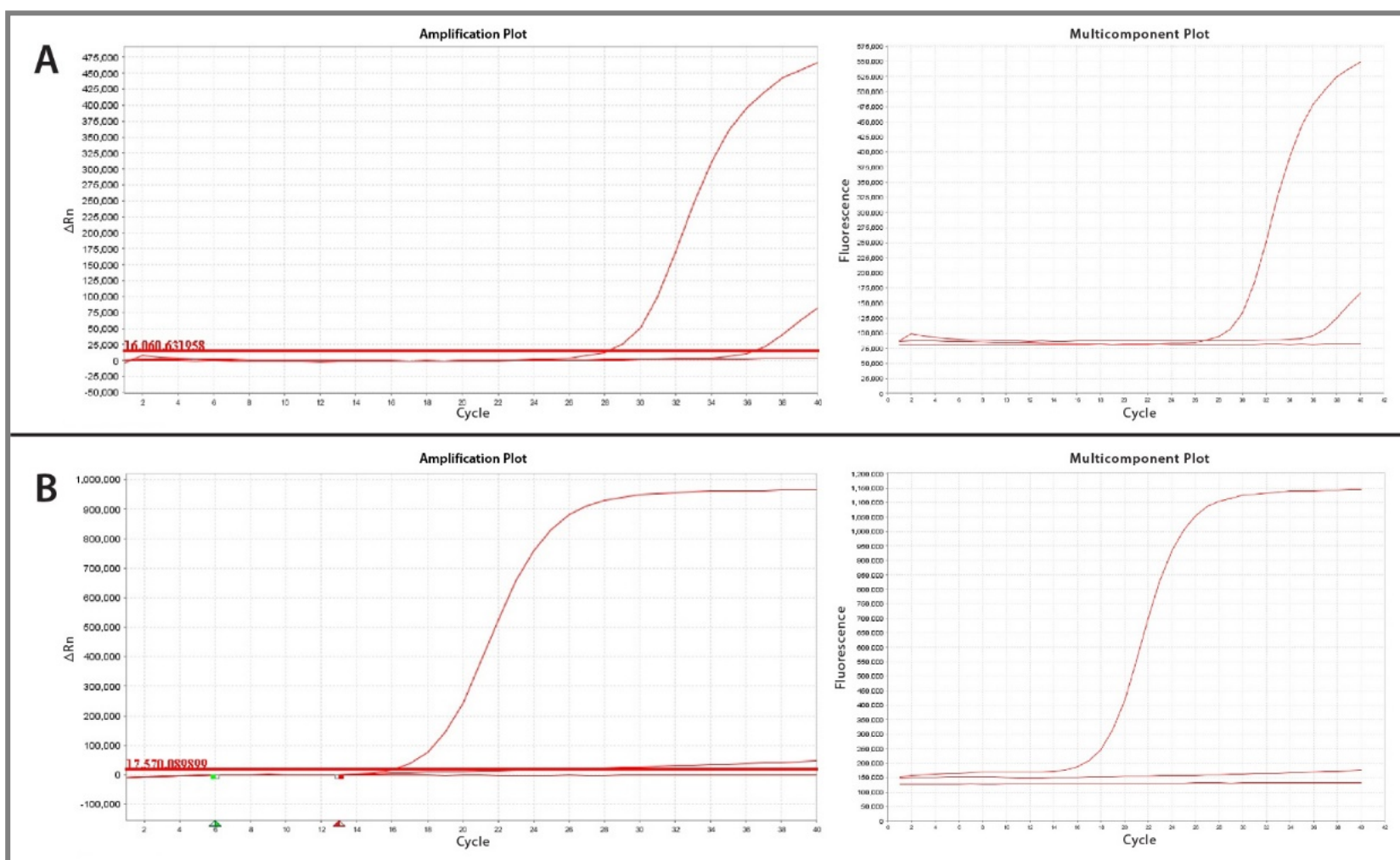Figure 13: Comparison between Real Amplification Signal (A) and Incorrect Baseline Setting (B)

Always check the signal displays in the Multicomponent Plot and the baseline is correctly set before concluding that an amplification trace is exponential. The Multicomponent Plot displays the complete spectral contribution of each dye and helps to review reporter dye signal for spikes, dips and/or sudden changes. Contact the equipment manufacturer or Fast Track Diagnostics for advice on how to correctly set up the baseline.

Limitations

- The use of the FTD SARS-CoV-2 as an *in vitro* diagnostic under the FDA Emergency Use Authorization (EUA) is limited to laboratories that are certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, to perform high complexity tests.
- FTD SARS-CoV-2 can be used only with the specimens listed in the Intended Use statement. Other specimen types have not been evaluated with this assay.
- The sample storage and shipment instructions provided are recommendations. Verification and validation data for specimen collection and handling (transport and storage) are not available.
- Remel M4RT® transport medium is not recommended for use with FTD SARS-CoV-2.
- This kit is a qualitative kit that does not provide a quantitative value for the detected pathogens in the specimen. There is no correlation between the Ct values obtained and the amount of pathogens in the specimen collected.
- Other parameters can lead to false positive, negative or invalid results related to patient conditions (use of antiviral therapy, patient age, patient history of respiratory infections, presence of symptoms and the stage of infection).
- Use of this kit should be limited to personnel trained in the technique of RT-PCR and in the use of FTD kits.
- The performance of the kit has been verified and validated using the procedures provided in the instructions for use only. Modifications to these procedures may alter the performance of the test.
- The performance of this kit has been evaluated for use with human specimen material only.
- This test shall not be the only element consulted for diagnosis or treatment decision. A specimen not detected cannot be presumed to be negative for this pathogen since results are dependent on several variables as explained above.
- Reliable results of this test require appropriate specimen collection as well as appropriate specimen and kit transport and storage and processing procedures. Failure to follow these procedures will produce incorrect results, leading to false positive and negative values or invalid results.
- Low levels of viruses can be detected below the limit of detection, but results may not be reproducible.
- Mutations within the regions of the targets for the virus detected by the kit may occur. As a consequence, primer and probe combinations may fail to detect the presence of this virus.

Conditions of Authorization for the Laboratory

The Fast Track Diagnostics™ SARS-CoV-2 test Letter of Authorization, along with the authorized Fact Sheet for Healthcare Providers, the authorized Fact Sheet for Patients and authorized labeling are available on the FDA website:

<https://www.fda.gov/medical-devices/emergency-situations-medical-devices/emergency-use-authorizations#covid19ivd>

To assist clinical laboratories using the FTD SARS-CoV-2 test, the relevant Conditions of Authorization are listed below:

- A. Authorized laboratories^[1] using the FTD SARS-CoV-2 test will include with test result reports, all authorized Fact Sheets. Under exigent circumstances, other appropriate methods for disseminating these Fact Sheets may be used, which may include mass media.
- B. Authorized laboratories using the FTD SARS-CoV-2 test will use the FTD SARS-CoV-2 test as outlined in the FTD SARS-CoV-2 Instructions for Use. Deviations from the authorized procedures, including the authorized instruments, authorized extraction methods, authorized clinical specimen types, authorized control materials, authorized other ancillary reagents and authorized materials required to perform the FTD SARS-CoV-2 test are not permitted.
- C. Authorized laboratories that receive the FTD SARS-CoV-2 test must notify the relevant public health authorities of their intent to run the test prior to initiating testing.
- D. Authorized laboratories using the FTD SARS-CoV-2 test will have a process in place for reporting test results to healthcare providers and relevant public health authorities, as appropriate.
- E. Authorized laboratories will collect information on the performance of the test and report to DMD/OHT7-OIR/OPEQ/CDRH (via email: CDRH-EUA-Reporting@fda.hhs.gov) and FTD Technical Support (tel: 1-877-229-8601 / email: support-ftd.team@siemens-healthineers.com) any suspected occurrence of false positive or false negative results and significant deviations from the established performance characteristics of the test of which they become aware.
- F. All laboratory personnel using the test must be appropriately trained in RT-PCR techniques and use appropriate laboratory and personal protective equipment when handling this kit and use the test in accordance with the authorized labeling.
- G. FTD, its authorized distributor(s) and authorized laboratories using the FTD SARS-CoV-2 test will ensure that any records associated with this EUA are maintained until otherwise notified by FDA. Such records will be made available to FDA for inspection upon request.

[1] For ease of reference, this letter will refer to, "Laboratories certified under the Clinical Laboratory Improvement Amendments of 1988 (CLIA), 42 U.S.C. §263a, to perform high complexity tests" as "authorized laboratories."

Performance Characteristics

Performance characteristics show the analytical and clinical performance data of FTD SARS-CoV-2. The analytical performance (analytical sensitivity, inclusivity and analytical specificity) was evaluated using the NucliSENS® easyMAG® (bioMérieux) extraction system with the Applied Biosystems® 7500 (ThermoFisher Scientific) PCR System and the VERSANT® kPCR Molecular System (Siemens Healthineers).

Limit of Detection – Analytical Sensitivity

Limit of detection (LoD) studies determine the lowest detectable concentration of SARS-CoV-2 at which greater or equal to 95% of all (true positive) replicates test positive.

To determine the LoD, a cultured virus (quantified) of an isolate from a US patient (USA-WA1/2020, vendor Zeptomatrix, catalog number 0810587CFHI) was serially diluted in simulated respiratory matrix (SRM).

A preliminary LoD was determined testing replicates of 2-fold serial dilutions of quantified cultured SARS-CoV-2 virus. The preliminary LoD data were used for the Probit (PROBability unitTs) Regression Analysis. Based on results provided by Probit analysis, an LoD confirmation study was performed. The LoD concentration from the Probit analysis (0.0023 TCID₅₀/mL) was confirmed using the NucliSENS® easyMAG® extraction system / Applied Biosystems® 7500 real-time thermocycler and the VERSANT® kPCR Molecular System.

NucliSENS® easyMAG® / Applied Biosystems® 7500

Twenty-three separate spiked samples were prepared by spiking the LoD concentration of the virus culture into SRM (in ESwab™ Liquid Amies medium [COPAN]). Each sample was extracted once (200 µL sample input, 55 µL elution volume) and tested in singlicate with 3 different lots of PP mix. The results of the LoD confirmation study are summarized in Table 8.

Table 8: Results of LoD Confirmation (NucliSENS® easyMAG® / Applied Biosystems® 7500)

PP Mix Lot	LoD (TCID ₅₀ /mL)	Tested	Total Detected	Detection Rate (%)	Mean Ct	SD
1	0.0023	23	22	95.7	38.0	0.9
2	0.0023	23	23	100.0	37.4	0.9
3	0.0023	23	22	95.7	37.4	0.8

Legend: TCID₅₀/mL = Median tissue culture infectious dose per milliliter, PP = Primer and probe, Ct = Cycle threshold, SD = Standard deviation

The LoD at 0.0023 TCID₅₀/mL of FTD SARS-CoV-2 in simulated respiratory matrix, as determined by the Probit analysis, was experimentally confirmed with an overall detection rate of 97.85%.

VERSANT® kPCR Molecular System

Twenty separate spiked samples were prepared by spiking the LoD concentration of the virus culture into SRM. Samples were prepared in SRM in three different types of transport media, including Universal Viral Transport (UVT, Becton Dickinson), Viral Transport Medium (VTM) prepared in-house according to CDC Protocol DSR-052-01, and ESwab™ Liquid Amies medium (COPAN). Samples were diluted to the LoD at 0.0023 TCID₅₀/mL as determined by the Probit Regression Analysis. Each sample was pretreated in a 1:2 dilution (350 µL sample with 350 µL VERSANT® Sample Preparation 1.0 Lysis Buffer). Samples were then extracted using the DAP FTD1 (SP1.0 500-15-10) protocol (500 µL input volume, 100 µL elution volume) and tested in singlicate. The results of the LoD confirmation study are summarized in Table 9.

Table 9: Results of LoD Confirmation (VERSANT® kPCR Molecular System)

Medium	LoD (TCID ₅₀ /mL)	Tested	Total Detected	Detection Rate (%)	Mean Ct	SD
SRM in UVT	0.0023	20	19	95.0	36.4	0.7
SRM in VTM	0.0023	20	20	100.0	36.5	0.9
SRM in ESwab™	0.0023	20	19	95.0	37.1	1.0

Legend: TCID₅₀/mL = Median tissue culture infectious dose per milliliter, PP = Primer and probe, Ct = Cycle threshold, SD = Standard deviation

The LoD at 0.0023 TCID₅₀/mL of FTD SARS-CoV-2 in simulated respiratory matrix in three different transport media was experimentally confirmed with an overall detection rate of greater than or equal to 95%.

Inclusivity – Analytical Sensitivity

The inclusivity *in silico* analysis was performed on all sequences available on the National Center for Biotechnology Information (NCBI) GenBank and Global Initiative on Sharing All Influenza Data (GISAID) databases. 1048 sequences (114 from GenBank and 934 from GISAID – sequences downloaded on 19 March 2020) were aligned against FTD SARS-CoV-2 primers and probes.

After excluding nonrelevant sequences, 901 sequences were determined to be appropriate for further analysis. *In silico* analysis concluded that FTD SARS-CoV-2 will detect all analyzed SARS-CoV-2 sequences in the NCBI GenBank (n=96) and in GISAID (n=805) databases.

FTD SARS-CoV-2 N gene assay detected all sequences from the GenBank database with a maximum of 1 mismatch (3 sequences) and 100% of all sequences obtained from the GISAID database with a maximum of 1 mismatch (11 sequences). The SARS-CoV-2 ORF1ab assay detected all sequences from the GenBank database without any mismatches and 100% of all sequences obtained from the GISAID database with a maximum of 1 mismatch (8 sequences). None of these mismatches were located at a critical position that would cause detection issues and are not predicted to impact assay performance.

Cross-Reactivity – Analytical Specificity

In Silico Analysis

A total of 39 bacterial/viral/fungal strains have been analyzed *in silico*. NCBI BLAST tool was used to check for cross-reactivity of the different primers and probes of the SARS-CoV-2 assay against the non-redundant nucleotide database. BLAST tool search default parameters were used except for the “organism.” The search was limited to using the taxonomy ID (taxid/txid) of the respective pathogen. Each primer and probe were compared against all available genome sequences of a certain taxid.

The results showed that in a few cases, the analyzed pathogens have more than 80% homology with any of the primers/probes designed for FTD SARS-CoV-2. Among those are SARS-CoV, *Legionella sainthelensi*, *Legionella spiritensis*, *Staphylococcus epidermidis*, and *Streptococcus salivarius*. These pathogens have been analyzed in more detail in a sequence alignment. No potential unintended cross-reactivity is expected based on this *in silico* analysis. Results are shown in Table 10.

Table 10: *In Silico* Cross-Reactivity Results^[a]

Pathogen	Taxonomy ID Included in BLAST Search	<i>In Silico</i> Analysis Result N gene Primer/Probe	<i>In Silico</i> Analysis Result ORF1ab Primer/Probe
<i>Influenza C virus</i>	NCBI:txid11552	< 80%	< 80%
<i>Parechovirus</i>	NCBI:txid138954	< 80%	< 80%
<i>Candida albicans</i>	NCBI:txid5476	< 80%	< 80%
<i>Corynebacterium diphtheriae</i>	NCBI:txid1717	< 80%	< 80%
<i>Legionella non-pneumophila</i> ^[b]	NCBI:txid445 exclude NCBI:txid446	83%	< 80%
<i>Bacillus anthracosis</i> (Anthrax)	NCBI:txid1392	< 80%	< 80%

Table 10: *In Silico* Cross-Reactivity Results^[a] (Continued)

Pathogen	Taxonomy ID Included in BLAST Search	<i>In Silico</i> Analysis Result N gene Primer/Probe	<i>In Silico</i> Analysis Result ORF1ab Primer/Probe
<i>Moraxella catarrhalis</i>	NCBI:txid480	< 80%	< 80%
<i>Neisseria elongata</i>	NCBI:txid495	< 80%	< 80%
<i>Neisseria meningitidis</i>	NCBI:txid487	< 80%	< 80%
<i>Pseudomonas aeruginosa</i>	NCBI:txid287	< 80%	< 80%
<i>Staphylococcus epidermidis</i>	NCBI:txid1282	83%	< 80%
<i>Streptococcus salivarius</i>	NCBI:txid1304	< 80%	91%
<i>Leptospira</i>	NCBI:txid171	< 80%	< 80%
<i>Chlamydia psittaci</i>	NCBI:txid83554	< 80%	< 80%
<i>Coxiella burnetii</i> (Q-Fever)	NCBI:txid777	< 80%	< 80%
<i>Staphylococcus aureus</i>	NCBI:txid1280	< 80%	< 80%
<i>Streptococcus pyogenes</i>	NCBI:txid1314	< 80%	< 80%
<i>Pneumocystis jirovecii</i>	NCBI:txid42068	< 80%	< 80%
Human coronavirus 229E	NCBI:txid11137	< 80%	< 80%
Human coronavirus OC43	NCBI:txid31631	< 80%	< 80%
Human coronavirus HKU1	NCBI:txid290028	< 80%	< 80%
Human coronavirus NL63	NCBI:txid277944	< 80%	< 80%
SARS-coronavirus	NCBI:txid694009 exclude HCoV-19 (taxid:2697049)	90%	95%
MERS-coronavirus	NCBI:txid1335626	< 80%	< 80%
Adenovirus (e.g., C1 Ad. 71)	NCBI:txid10509	< 80%	< 80%
Human Metapneumovirus (hMPV)	NCBI:txid162145	< 80%	< 80%
Human Parainfluenza virus 1–4	NCBI:txid12730, NCBI:txid11216, NCBI:txid1979160, NCBI:txid1979161	< 80%	< 80%
Influenza A virus	NCBI:txid11320	< 80%	< 80%
Influenza B virus	NCBI:txid11520	< 80%	< 80%
Enterovirus A–D (e.g., EV68)	NCBI:txid138948, NCBI:txid138949, NCBI:txid138950, NCBI:txid138951	< 80%	< 80%
Human Respiratory syncytial virus	NCBI:txid11250	< 80%	< 80%
Rhinovirus	NCBI:txid147711, NCBI:txid147712, NCBI:txid463676	83%	< 80%
<i>Chlamydia pneumoniae</i>	NCBI:txid83558	< 80%	< 80%

Table 10: *In Silico* Cross-Reactivity Results^[a] (Continued)

Pathogen	Taxonomy ID Included in BLAST Search	<i>In Silico</i> Analysis Result N gene Primer/Probe	<i>In Silico</i> Analysis Result ORF1ab Primer/Probe
<i>Haemophilus influenzae</i>	NCBI:txid727	< 80%	< 80%
<i>Legionella pneumophila</i>	NCBI:txid446	< 80%	< 80%
<i>Mycobacterium tuberculosis</i>	NCBI:txid1773	< 80%	< 80%
<i>Streptococcus pneumoniae</i>	NCBI:txid1313	< 80%	< 80%
<i>Bordetella pertussis</i>	NCBI:txid520	< 80%	< 80%
<i>Mycoplasma pneumoniae</i>	NCBI:txid2104	< 80%	< 80%

[a] Results denote the percent coverage given by BLAST analysis for all primers and probes, in case homology was greater than 80%, it shows the highest value.

[b] Among the *Legionella non-pneumophila* strains, *Legionella sainthelensi* and *Legionella spiritensis* revealed a homology greater than 80% for given primers and probe.

In addition, the concerned organisms have been tested *in vitro* for potential cross-reactivity and no amplification was observed (see *In Vitro* Analysis section).

***In Vitro* Analysis**

A total of 32 bacterial/viral RNA/DNA cultures or samples have been tested *in vitro* for cross-reactivity. A total of 5 pools were generated by spiking a maximum of 5 organisms into Eswab™ Liquid Amies medium (COPAN). In addition, a pool of 5 nasal fluid patient samples was generated. Each pool was extracted in triplicate and tested with FTD SARS-CoV-2. Six cultures have been extracted in triplicates and tested separately without pooling. A culture was not available for human coronavirus HKU1, but a patient sample was extracted and tested. In addition, 6 genomic previously-acquired RNA/DNA samples were tested for PCR.

Table 11 lists all tested organisms and the respective tested concentration. For some organisms the concentration was unknown, they were tested previously by RT-PCR using FTD Respiratory pathogens 21 (CE-IVD) and Ct values are indicated. For others, the Ct values given by the vendor are indicated.

Table 11: List of Pools or Individual Pathogens Tested *In Vitro*

Pool ID / Sample ID	Interfering Organism	Tested Concentration	Unit	Result
P1	Adenovirus 71	1.13E+06	TCID ₅₀ /mL	No cross-reactivity
	Human Parainfluenza virus 2	1.13E+06	TCID ₅₀ /mL	No cross-reactivity
	Enterovirus	2.00E+06	TCID ₅₀ /mL	No cross-reactivity
	Rhinovirus	1.13E+07	TCID ₅₀ /mL	No cross-reactivity
	<i>Streptococcus pneumoniae</i>	1.24E+08	cop/mL	No cross-reactivity

Table 11: List of Pools or Individual Pathogens Tested *In Vitro* (Continued)

Pool ID / Sample ID	Interfering Organism	Tested Concentration	Unit	Result
P2	Human Metapneumovirus (hMPV) A	6.82E+05	TCID ₅₀ /mL	No cross-reactivity
	Parainfluenza virus 3	1.13E+06	TCID ₅₀ /mL	No cross-reactivity
	Parainfluenza virus 4	2.00E+05	TCID ₅₀ /mL	No cross-reactivity
	Human Metapneumovirus (hMPV) B	1.56E+05	TCID ₅₀ /mL	No cross-reactivity
	Parainfluenza virus 1	2.00E+05	TCID ₅₀ /mL	No cross-reactivity
P3	Human coronavirus OC43	1.41E+04	TCID ₅₀ /mL	No cross-reactivity
	Human coronavirus NL63	8.39E+03	TCID ₅₀ /mL	No cross-reactivity
	Influenza B virus	Unknown (Ct 20.4) ^[a]	N/A	No cross-reactivity
P4	<i>Chlamydomphila pneumoniae</i>	2.01E+04	cop/mL	No cross-reactivity
	<i>Haemophilus influenzae</i>	5.71E+02	CFU/mL	No cross-reactivity
	<i>Legionella pneumophila</i>	1.00E+03	CFU/mL	No cross-reactivity
	<i>Bordetella pertussis</i>	1.43E+03	CFU/mL	No cross-reactivity
P5	<i>Mycoplasma pneumoniae</i>	Unknown (Ct 21) ^[a]	N/A	No cross-reactivity
	Human Respiratory syncytial virus (HRSV-A)	3.14E+06	PFU/mL	No cross-reactivity
	Human Respiratory syncytial virus (HRSV-B)	3.26E+05	TCID ₅₀ /mL	No cross-reactivity
P6	Nasal Fluid Pool	N/A	N/A	No cross-reactivity
S7	<i>Pneumocystis jirovecii</i>	1.64E+08	cop/mL	No cross-reactivity
S8	<i>Streptococcus salivarius</i>	5.91E+07	cop/mL	No cross-reactivity
S9	<i>Staphylococcus epidermidis</i>	8.05E+06	cop/mL	No cross-reactivity
S10	<i>Legionella sainthelensi</i>	2.29E+08	cop/mL	No cross-reactivity
S11	<i>Legionella spiritensis</i>	4.08E+08	cop/mL	No cross-reactivity

Table 11: List of Pools or Individual Pathogens Tested *In Vitro* (Continued)

Pool ID / Sample ID	Interfering Organism	Tested Concentration	Unit	Result
S12	<i>Streptococcus pyogenes</i>	6.82E+08	cop/mL	No cross-reactivity
S13	<i>Mycobacterium tuberculosis</i>	1.00E+06	cop/mL	No cross-reactivity
S14	Influenza A virus	1.00E+06	cop/mL	No cross-reactivity
S15	SARS-coronavirus	Unknown (Ct 18) ^[a]	N/A	No cross-reactivity
S16	MERS-coronavirus	Unknown (Ct 29) ^[a]	N/A	No cross-reactivity
S17	Human coronavirus 229E	Unknown (Ct 28) ^[a]	N/A	No cross-reactivity
S18	Human coronavirus HKU1	Unknown (Ct 14) ^[a]	N/A	No cross-reactivity

[a] Ct value given by supplier or tested with FTD Respiratory pathogens 21 (CE-IVD).

Legend: TCID₅₀ = Median tissue culture infectious dose, cop/mL = copies per milliliter, CFU/mL = Colony-forming unit per milliliter, PFU/mL = Plaque-forming unit per milliliter

Clinical Performance

The performance of FTD SARS-CoV-2 was evaluated using the NucliSENS® easyMAG® extraction system / Applied Biosystems® 7500 real-time thermocycler and the VERSANT® kPCR Molecular System.

NucliSENS® easyMAG® / Applied Biosystems® 7500

The performance of FTD SARS-CoV-2 was established using 80 nasopharyngeal swabs, collected from male and female adult patients with signs and symptoms of an upper respiratory infection.

A total of 44 positive specimens and 36 negative specimens were tested with FTD SARS-CoV-2 using the NucliSENS® easyMAG® extraction method and the Applied Biosystems® 7500 Real-Time PCR System. The clinical performance study was evaluated by comparing FTD SARS-CoV-2 results to the FDA EUA-authorized RT-PCR test. During the analyses, 6 negative samples were excluded; 1 was excluded because of a protocol deviation and 5 did not pass criteria for a valid run (IC failure). Results are displayed in Table 12.

Table 12: Positive and Negative Percent Agreements (NucliSENS® easyMAG® / Applied Biosystems® 7500)

Agreements between FTD SARS-CoV-2 and FDA EUA RT-PCR using Nasopharyngeal Swabs (n=74)		FDA EUA RT PCR Test	
		Positive	Negative
FTD SARS-CoV-2	Positive	44	0
	Negative	0	30
Positive Percent Agreement	100% (44/44) (95% Confidence Interval: 91.97, 100)		
Negative Percent Agreement	100% (30/30) (95% Confidence Interval: 88.65, 100)		

The results showed 100% positive percent agreement (95% Confidence Interval: 91.97–100) and 100% negative percent agreement (95% Confidence Interval: 88.65–100) between FTD SARS-CoV-2 and the FDA-authorized RT-PCR test, for the detection of SARS-CoV-2 in nasopharyngeal swabs.

VERSANT® kPCR Molecular System

The performance of FTD SARS-CoV-2 was determined using 109 nasopharyngeal swabs, collected from male and female patients with signs and symptoms of an upper respiratory infection.

A total of 60 positive specimens and 49 negative specimens, as determined using an FDA EUA-authorized RT-PCR test, were tested with FTD SARS-CoV-2 using the VERSANT® kPCR Molecular System. The clinical performance study was evaluated by comparing FTD SARS-CoV-2 results to the FDA EUA RT-PCR test. Each sample was pretreated in a 1:2 dilution (350 µL sample with 350 µL VERSANT® Sample Preparation 1.0 Lysis Buffer) prior to extraction. Results are displayed in Table 13.

Table 13: Positive and Negative Percent Agreements (VERSANT)

Agreements between FTD SARS-CoV-2 and FDA EUA RT-PCR test using VERSANT® kPCR Molecular System (n=109)		FDA EUA RT PCR Test	
		Positive	Negative
FTD SARS-CoV-2	Positive	60	3*
	Negative	0	46
Positive Percent Agreement	100% (60/60) (95% Confidence Interval: 93.98, 100)		
Negative Percent Agreement	93.88% (46/49) (95% Confidence Interval: 83.48, 97.90)		

*2 of 3 discordant samples (positive by FTD SARS-CoV-2) were positive by sequencing.

The results showed 100% positive percent agreement (95% Confidence Interval: 93.98–100) and 93.88% negative percent agreement (95% Confidence Interval: 83.48–97.90) between FTD SARS-CoV-2 and the FDA EUA-authorized RT-PCR tests, for the detection of SARS-CoV-2 in nasopharyngeal swabs. Two of the three discordant results (positive with FTD SARS-CoV-2) were also positive using sequencing methods.

Troubleshooting

Table 14 describes a non-exhaustive list of control errors that a user may observe with FTD SARS-CoV-2 and suggested corrective actions.

Table 14: Control Errors

Observation	Possible Cause	Corrective Action
Positive control or external RNA positive control does not amplify	Incorrect programming of the thermocycler temperature profile.	Compare temperature profile to IFU.
	Incorrect configuration of the PCR run.	<ul style="list-style-type: none"> Confirm reagents were added in the correct sequence; repeat the PCR, if necessary. Check calibration of pipettes.
	Incorrect handling of the positive controls.	Inadequate or no vortexing, or control was not adequately thawed at room temperature.
	Storage conditions for one or more product components did not comply with the instructions or the FTD kit has expired.	Check storage conditions and expiration date on the kit box. Discard the kit if necessary.
Weak or no signal of the internal control	PCR conditions do not comply with protocol.	Ensure extraction and amplification workflow was performed as described. Repeat analysis, if necessary. In case the problem persists, consider the presence of interfering material in your samples.
	Amplification of IC was inhibited or the extraction of the IC was inadequate.	
Amplification in the negative control	Contamination during PCR plate set up or during extraction.	<ul style="list-style-type: none"> Repeat PCR plate set up with new reagents, samples and controls. Repeat extraction procedure with new reagents. To avoid contamination from the PC, pipette the positive control last. Decontaminate the workspace and instruments after each use.

If the problem persists, note the error and contact technical support, go to www.fast-trackdiagnostics.com.

Technical Assistance















For customer support, contact the local technical support provider or distributor or refer to the Technical Support section of the Fast Track Diagnostics website at www.fast-trackdiagnostics.com.

References

1. The Institut Pasteur website:
<https://www.pasteur.fr/en/medical-center/disease-sheets/covid-19-disease-novel-coronavirus#symptoms>. Accessed March 2020.
2. Center for Disease Control (CDC) website:
<https://www.cdc.gov/urdo/downloads/SpecCollectionGuidelines.pdf>. Accessed March 2020.
3. CDC Interim Guidelines for Collecting, Handling, and Testing Clinical Specimens from Patients Under Investigation (PUIs) for 2019 Novel Coronavirus.
<https://www.cdc.gov/coronavirus/2019-nCoV/guidelines-clinical-specimens.html>. Accessed May 2020.
4. World Health Organization (WHO). Laboratory testing for coronavirus disease 2019 (COVID-19) in suspected human cases: Interim guidance, 2 March 2020.
5. Siemens Healthineers. SMN 11222579, *VERSANT® kPCR Molecular System Sample Preparation (SP) Application Guide*.
6. Siemens Healthineers. SMN 11223551, *VERSANT® kPCR AD Application Guide (Thermo Fisher Scientific QuantStudio 5 Dx)*.

Definition of Symbols

This section describes all symbols used to convey product labeling description, use or handling information on components or unit-of-sale packaging.

Symbol	Definition	Symbol	Definition
	<i>In vitro</i> diagnostic medical device		Contains sufficient for <n> tests
	Catalog number		Batch code
	Manufacturer		Use-by date
	Date of manufacture		Keep away from sunlight
RxOnly	Prescription device (US Only) Applies only to United States IVD assays. CAUTION: Federal (USA) law restricts this device to sale by or on the order of a licensed healthcare professional.	YYYY-MM-DD	Date format (Year-Month-Day)
	Consult instructions for use	YYYY-MM	Date format (Year-Month)
	Caution/Warning		Store upright
	Temperature limit		Irritant
			Made in Luxembourg

Legal Information

FTD and all associated marks are trademarks of Fast Track Diagnostics Luxembourg S.à r.l. or its affiliates. All other trademarks and brands are the property of their respective owners.

© 2020 Fast Track Diagnostics. All rights reserved.

Siemens Healthineers Headquarters

Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
siemens-healthineers.com



Legal Manufacturer

Fast Track Diagnostics Luxembourg S.à r.l.
29, rue Henri Koch
L-4354 Esch-sur-Alzette
Phone: Luxembourg +011-352-281098-217 / USA 1-877-229-8601
Email: support-ftd.team@siemens-healthineers.com